EASY BAND

USER MANUAL

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


MANUAL DE UTILIZACION

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ELECTROSTIMULATOR

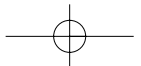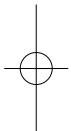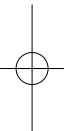SUMMARY

ENGLISH

English

INDEX

-Accessories	2
-Precautionary measures	3
-Technical features	5
-Electrode and battery removal	
-Description of the Band device	6
-Description of the Band Display	7
-How to use the Band	8
-How to remove the stimulator from the belt	9
-Battery replacement	
-Electrode and battery removal	
-How to put on the Band	
-How to place the stimulator in the Band	
-First use	
-What to do in case the current feels unpleasant	
-Useful advises	
-How to pause the program	
-Automatic functions	
-Automatic switching-off	
-Electrode contact	
-Empty battery	
-Intensity of stimulation	
-Special functions	
-Keyboard lock key	
-Common questions	11
-Guarantee terms	12

ACCESSORIES

Dotazione

- Electro stimulator
- Carry case
- Abdominis belt
- 3 Battery stylos 1,5 V AAA
- Manual

Accessories on request

- Gel

This stimulator is supplied complete of cables and electrodes. When you open the packing verify that the simulator equipment is complete. If some of the above mentioned accessories are missing, you should immediately contact the shopkeeper. Also check the stimulator and electrodes integrity.

Electrodes: The electrodes are usable by only one user they are self-adhesive and do not need gel. The electrode cable is female. They are branded CE complied with MDD 93/42/EEC requirements.

All the information can be subject to changes without notice.

The symbols mean that the apparatus is in accordance with medical CE requirements and type BF.

PRECAUTIONARY MEASURES

This symbol warns you to read the operating instructions prior to use

Contraindications

- Do not stimulate the neck region (carotid sinus).
- Patients who have a demand type cardiac pacemaker.
- Stimulation should be not used on cancer patients (previous permission received by the oncologist).
- Do not stimulate the cerebral zone.
- Do not use on undiagnosed pains until aetiology is established.
- Do not use on sores and dermatological diseases.
- Do not use on acute traumas.
- Do not stimulate recent cicatrix.
- Do not use during pregnancy.
- During driving car.
- It is absolutely forbidden to stimulate the ocular and genital zone
- It is also forbidden the use of the stimulator in the water (it is not protected against the water seepage). Not to use the stimulator in environment where are present combustible gas, inflammables, or oxygen.
- It is also required use with caution in case of capillary fragility because an excessive stimulation could cause the break of more quantity of capillaries. Anyway it is always preferable to contact your doctor before using.

Adverse possible effects

Skin irritation can occur in rare case. In case of allergic reactions by the use of gel, stop the treatment and contact the specialist doctor.
If during the stimulation treatment tachycardia or extrasistolia should appear, it is better to stop the treatment the sensation will pass immediately.

Precautions for using

- Do not make any modifications or repairs by yourself. Authorized centres can only execute the repairs.
- Do not give a shock to the product by clashing, overturning or shaking. It may cause a trouble or accident in a long run although no immediate trouble occurs.
- Store the electro stimulator at a place free from a splash or water. Avoid installing at places subject to undue humidity, temperature, dust or places exposed to direct sunbeam.
- Do not use thinner, gasoline, kerosene, cleaning powder or chemicals; for cleaning the machine it is sufficient a soft duster
- Take away the battery from the electro stimulator in case of prolonged no use.
- After a long period of no use it is advisable to check the correct working.
- For cleaning the machine and the cables, it is sufficient a soft duster lightly humid.

Warnings

- The stimulator must be used only for transcutaneous neuro-stimulation.
- The neurostimulation apparatus should be used only under the continued supervision of a physician.
- The stimulator should be only used with skin electrodes intended for nerve and muscle stimulation.
- The user must keep the device out of reach of children.
- Electronic monitoring equipment (such as ECG monitors and ECG alarms) may not operate properly when stimulation is in use
- Stimulus delivered by this device may be sufficient to cause electrocution: Electrical current of this magnitude must not flow through the thorax because it may cause a cardiac arrhythmia.
- Do not use the apparatus in case of any health problem without consulting your doctor.
- In case of damage of the covering apparatus, this has to be replaced.
- It is advisable not to use the apparatus together with other electronic equipment above all those used for support of vital functions.
- The device must be used only by person with mental competence.
- Do not use the stimulator near short wave devices.

Use destination

Read the operating manual carefully in order to avoid to carry out installation, use and maintenance non adequate that could damage the apparatus the operator and the user.

The stimulators are used in the following operating environment and field:

- home;
- surgeries hospital;
- physiotherapeutic surgery;
- rehabilitation centres;
- treatments of pain in general
- in any type of sequential electro stimulation,
- Aesthetic and sports use

The use of the stimulator is allowed to patients (rightly informed about the methods), doctors and paramedic. The presence of clinic stimulation, makes this product a medical device. The certification cover the clinic applications.

TECHNICAL FEATURES

Number of channels:	2.
Out current:	90mA
Wave form:	Square, biphasic, symmetric, compensated impulse.
Running frequency:	1-110Hz.
Impulse amplitude:	50-300s.
Supply:	3 batteries of 1,5V AAA type LR03.
Conditions of conservation and transport:	
Temperature:	10-40°C.
Air humidity:	20-65%.
Air pressure:	700hPa-1060hPa.
Weight:	95 gr.
Classification:	Inner supply type BF.
Certification:	Certificate CE MDD (Medical Device Directive).

The name of producer is reported in the last page

Device removal

Don't throw the battery or electrodes out but recycle or get rid of them appropriately following the removal rules in force in your country.


DESCRIPTION OF THE BAND DEVICE

ON / OFF key

- Press the key to switch-on the unit, press for 3 seconds to switch-off.
- Be switch also used during the treatment as a stop button

Program & LOCK key

- Press the key to select the following program.
- Press for 3 seconds to lock the keys except for the ON / OFF key.
- To unlock the keys, press for 3 seconds.

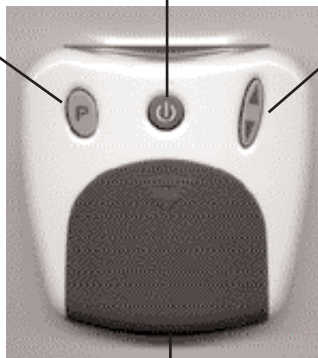Intensity key

- Press the key with the arrow ▲ to increase the out intensity
- Press the key with the arrow ▼ to decrease the out intensity.

Battery compartment


- Push down the cover and put 3 batteries of 1,5V AAA.
- Respect the indicated polarity.

Display

It gives different indications during the stimulation. The device has a big display that continuously gives useful information about the parameters in use during the stimulation.

DISPLAY

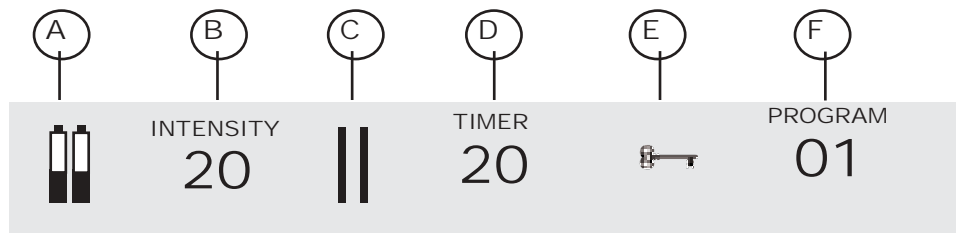English

The device has a big display, which continuously gives information about the different parameters in use.

A Empty battery indicator:

The indicator lightens when the batteries are almost empty.

B Intensity numeral indicator:

Shows the used intensity level included from 0 to 20 scale.

C Pause indicator:

The symbol II lightens when the stimulator is in pause.

D Timer:

Shows the remaining time in minutes.

E Lock Key indicator:

Shows that the keyboard is locked

F Program indicator:

Shows the selected program.

HOW TO USE THE BAND

How to remove the stimulator from the belt

Shut off the stimulator.

Take out the stimulator from the belt by pushing out the top.

When the unit is disconnected from the magnetic contact lift it upwards.

Battery replacement

Shut off the stimulator. Remove the battery cover by sliding it downwards and take out the batteries. Place 3 alkaline batteries of 1,5V AAA type LR03 and close the cover. Make sure that the new batteries are correctly placed with respect to the polarity (note marks for + and-)

Don't throw the battery or electrodes out but recycle or get rid of them appropriately following the removal rules in force in your country.

How to put on the BAND

The Band has been created to be adapted to a large number of different sizes. The special material allows an optimal contact to the skin that guarantees an efficacious and comfortable stimulation.

Dampen the electrode and the skin uniformly with water or gel and place the side with the black electrodes directly on the skin.

It is advisable to tighten the band so that to assure a perfect contact with the skin and muscle.

How to place the stimulator in the Band

Shut off the stimulator. To place the stimulator in the band, simply insert it in the opposite compartment with the display upwards. A "click" warns you that the unit is well linked.

First use

1. Press the ON / OFF key to turn-on the stimulator. The display lights and shows the last performed program. The intensity indicator shows 0.
2. Select the program you wish to perform by pressing the P key.
3. Increase the used channel intensity pressing the key with the arrow ▲ . Increase the intensity until the wished level, always maintaining a comfortable feeling. To decrease the intensity press the arrow ▼ .
4. Each program has a predefined duration. The intensity automatically returns to zero when the program ends.
5. Every time you wish to turn-off the unit, press the ON / OFF key for 3 seconds.

What to do in case the current feels unpleasant

If you feel pinches instead of muscular contraction, stop the program and follow the above indications:

- Take off the Band and be sure that the electrode zone is uniformly wet.
- Check that the electrodes are not damaged.
- When you fasten the Band, be sure that the electrodes have a good contact with the skin and the Band is good dressed around the stomach.

Useful advises

- Be sure that the conducting electrodes are always wet before using.
- If skin irritation should occur, treatment should temporarily discontinued. Often this effect is because the wet electrodes dry.
- A decrease of electric impulse effect could occur because an excessive perspiration. Wipe your skin with a cleaned cloth.
- Take a pause of at least 5 hours between two working sessions.

How to pause the program

You can pause the program by pressing the ON / OFF switch. Start the stimulation by pressing the On / OFF switch again. The program continues where the stimulation was shut off.

Automatic functions

Automatic Switching-off

If the electric stimulator is not in use for more than 5 minutes and the intensity in both channels is on 0, the device automatically switches-off.

Electrode contact

If the electrodes are not correctly placed on the skin or the cables are not well connected, the intensity of the channels is automatically shut off to 0.

In case the electrodes are not correctly placed before starting a treatment, it is not possible to increase the channel intensity.

This is a security measure to preserve the user from unpleasant possible effects.

Empty battery

The battery indicator lights up on the display when the batteries are almost empty. As long as stimulation feels the same as usual you can continue the stimulation. It is time to replace the batteries when stimulation starts diminishing or the stimulator shut off.

Intensity of stimulation

The intensity value for obtaining a particular contraction is completely subjective. It can depend on a series of factors like as the electrode position, the adipose stratum, the perspiration etc. For these reasons the same current intensity value can give different sensations from a person to another from a day to another or furthermore you could feel more contraction on the left limb rather than the right one.

You have to regulate it on the bases of your sensations.

Special Functions

Keyboard lock key

Once reached the wished intensity level, it is possible to lock the keys except for the ON / OFF switch by pressing the LOCK key for 3 seconds. This function will avoid you their accidental use during the treatment. The keyboard lock is released when the key is pressed for three seconds again.

Note: You cannot change the intensity value during the rump up and ramp down in all intermittent programs that foresees cycles of work and recovery.

COMMON QUESTIONS

Why does the intensity return immediately to 0 when increasing its value?

The unit has a control of current supply. In case the circuit is open the intensity is automatically shut off. It is therefore necessary, before starting any treatment, to be sure that the cables are well connected and the electrodes completely placed on the skin.

What to do when the discharged battery indicator appears?

You can continue your stimulation until it starts diminishing; then turn off the unit, disconnect the electrode cables and change the batteries.

Is it possible to use the electric stimulator during the menstrual cycle?

Possible interferences, like anticipations, delays, accentuations or reductions of the cycle can be subjective and variable. It is however advisable to avoid treatments in the abdominal zone during the cycle.

Is it possible to use the electric stimulator during lactation?

Until now it has not found collateral effects in relation to the electric stimulator use during the lactation. It is however advisable not to treat the thoracic zone.

Does the presence of dermatological diseases like psoriasis, urticaria, etc. prevent the electric stimulator use?

Yes, it is advisable not to treat epidermic zones, which present important dermatological diseases.

How long after is it possible to see the first results?

The aesthetic results obtained are always subjective. It is however possible to state that you can achieve good results about tonicity after 15 days with 3 or 4 sessions a week; for lypolysis and drainage treatment it is necessary at least 40 days. Better and speedier results are obtained if treatments are combined with a good physical activity and a correct life style.

GUARANTEE TERMS

This product is warranted to the first user for a 24-month period from the purchase date against material defects or manufacturing, provided that it is used properly and maintained in the normal efficiency conditions.

To enjoy the warranty service, the user must observe the following warranty articles:

1-The products must be delivered for repairing by the User in their original packaging, the expenses are charged to the User.

2-The product warranty depends on the presentation of fiscal document (cash voucher, fiscal receipt or selling invoice) and certifying the purchase date of the product.

3-The repairing will not have effect on the original date of warranty expiration and will not mean the renewing or extension of the warranty.

4-When during the repairing will not be found any defect, will be in any case charged the expenses related to the checking time.

5-The warranty doesn't have effect in case of damages by bits, falls, improper or incorrect use, not original external supply, accidents, alterations, changes / removal of warranty seals and / or product infringement. The warranty doesn't cover the damages caused during the transportation due to the use of inadequate boxing (see art. 1).

6-The warranty doesn't cover the impossibility of using the product, other occasional or consequential expenses or other purchaser expenses.

ATTENTION: Before returning the machine to be repaired, it is recommended to carefully read the user manual instructions. In case of device breakdowns contact the shopkeeper before. The fabricant can make at any moment all the necessary modifications to improve the product aesthetic and quality.

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Realizado da:

DOMINO s.r.l.
Tel. +390438793052
Fax +390438796463
info@domino.tv.it

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**This product as been manufactured by:**

DOMINO s.r.l.
Tel. + 390438793052
Fax +390438796463
E-mail: info@domino.tv.it

This product has been manufactured according to the technical regulations in force and it has been subjected to certifications according to 93/42/CEE directive for medical devices in order to assure the product's security. The user must follow all prescriptions and restrictions shown in the operating manual in order to keep the maximum security level of the product. The manufacturer declines all responsibility for every different use from what indicated in this operating manual. All complete or partial reproductions of texts and/or pictures of this operating manuals made by any electronic or mechanic equipments are strictly forbidden without written permission of the manufacturer.

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Accordance with 72/23/CEE low voltage and 89/336/CEE, 93/42/CEE with the medical device directive requirement. The number of the notified authority is 0476

This symbol warns you to read the operating instruction prior to use

This symbol means that this device is a BF pattern

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