*Tecnovita*

by **BH**

Instrucciones de uso y garantía  
Operating instructions and guarantee  
Mode d'emploi et garantie  
Instruções de utilização e garantia  
Manuale d'uso e garanzia  
Bedienungsanleitung und Garantie

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**902 170 258**

Thank you for purchasing our e-shorts YR40.  
All of our products undergo rigorous testing to ensure maximum quality.  
For a correct use of the unit, first read the whole manual before using it.

## INTRODUCTION

**EMS, electric muscular stimulation** is a widely recognized method that it has been used on the sports field for a long time. The e-short is both easy and comfortable to use. This product is specially designed to help tone thighs and gluteus. It allows enjoying the exercise and increasing muscle capacity by developing strength and improving tone effortlessly. Thanks to its design, the e-short can tone the abdominal muscles anytime, anywhere.

## GENERAL SAFETY INSTRUCTIONS


1. Before using the appliance make sure it is in good working order. If you detect any faults get in touch with your local dealer.
2. Children or disabled people without suitable adult supervision must not use the unit.
3. Pregnant and in post partum period women must not use it.
4. Do not use if you suffer from heart disease or other type of disease.
5. Do not use the appliance if you have cancer.
6. Do not use if you are prone to attacks, for example, epilepsy.
7. Do not use the appliance after surgical operations in which increased muscle tension may affect the recovery process
8. Do not use directly on open wounds or varicose veins.
9. Do not use if you have any implanted electrical device, such as a pacemaker.
10. Do not use if you have any metal implants.
11. Do not use on broken or inflamed skin.
12. If you suffer from any serious illness, especially if you suspect that you may have high blood pressure or problems with blood clotting, consult your doctor before using the appliance.
13. Consult your doctor before using this product if you have any type of skin disease.
14. If you suffer from chronic pain of unknown origin, regardless of where on your body, consult your doctor before using the appliance.
15. People undergoing medical treatment should consult a doctor before using this appliance.
16. Diabetics should consult their doctor before using this appliance.

17. People with constant skin irritation resulting from prolonged stimulation by the same electrode should consult a doctor.
18. Consult your doctor if you have any problems resulting from the stimulation treatment.
19. Stop using the appliance if you experience any discomfort during use.
20. Do not use after consuming alcohol.
21. If the appliance malfunctions or if you feel discomfort or pain, discontinue use immediately.
22. Do not use while sleeping, driving or while using some other device.
23. Keep the appliance away from heat sources and do not use close (<1 m) to medium or microwave devices (e.g. mobile phones) because they can cause power surges.
24. Do not use when there are high frequency surgical devices connected simultaneously.
25. Do not allow the electrodes to meet metal objects, e.g. belts or collars, during operation.
26. Do not use near flammable, gaseous or explosive substances.
27. Switch the appliance off before removing it to prevent unnecessary irritation.
28. Never use the appliance without first applying plenty of electro-stimulation gel. If you do, then the stimulating effect will not be transmitted correctly and it can damage the appliance.
29. Do not apply enriching creams beforehand, as these cause the electrodes to wear out faster or may cause unwanted current surges.
30. Do not use the unit in damp/humid surroundings and prevent it from getting wet. Do not touch the device with wet hands while it is operating. Under no circumstances, this product is to be used in a bath or shower.
31. Do not expose the appliance to direct sunlight or high temperatures.
32. Do not allow the appliance to meet water or any other liquid. Do not immerse the appliance in liquid to clean it.
33. If you notice that the appliance is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

## DESCRIPTION

### REMOTE CONTROL

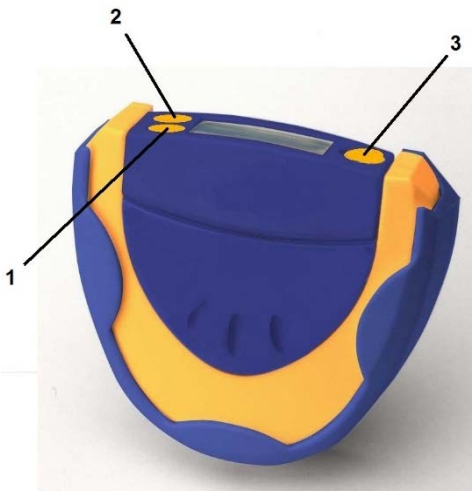- 1 – ON/Intensity + (10 levels)
- 2 – OFF/Intensity – (10 levels)
- 3 – Mode selection (6 modes)

MODE	DESCRIPTION	TRAINING SESSION
1	5 pulses per second	Aerobic fat burning
2	1 pulses per second	Aerobic fat burning
3	3 sec. ON 2 sec. OFF	Aerobic fat burning
4	Mixed sequence	Aerobic and toning
5	Advanced mixed sequence	Aerobic and toning
6	Continuously ON at low frequency	Toning

## OPERATING INSTRUCTIONS.

### BEFORE USE

1. Remove the control from its holder;
2. Open the cover sliding it downward;

3. Insert AAA type batteries (2 units) and close the cover;
4. Place the remote back in the remote holder.

## APPLICATION

1. Wear the pants and adjust the size through the fastening part (Velcro);
2. Put the device into the device holder on the pants with two snaps on buttons;
3. Push ON/HI button to turn on the device. To switch the device off press the OFF/LOW button.
4. Pressing the “MODE” button repeatedly to select the mode you wish. The LCD display will show the mode which you have selected. (This button will only be effective when the device is connected to the pants and in touch with the skin).

## INTENSITY ADJUST

Push the ON/HI or OFF/LO button to carefully increase or decrease the intensity. The LCD display will show the intensity level which you have selected. (This button will only be effective when the device is connected with the attachment and it touch on the skin) You will notice a tingling sensation prior to the muscle contracting. Once the muscle is contracting will be pleasant and relaxing. Adjust the intensity according to your comfort level. You should feel a strong but comfortable contraction of the muscles.

When the conductive area (SILICON PADS) of the pants does not touch the skin, the intensity can only be set up to “1”.

If the pants are removed during use, the intensity is automatically reset to “1”. If the attachment is not placed back into contact with the skin, the unit will turn off after 30 seconds. The unit has an auto-timer. It will shut off after 10 minutes of massage. If you change the massage mode within the 10 minutes, the timer will count-down from the original leftover time.

## POSSIBLE PROBLEMS

PROBLEM	CAUSE	SOLUTION
THE SCREEN DOES NOT SWITCH ON	<ul style="list-style-type: none"> <li>- THERE ARE NO BATTERIES INSTALLED</li> <li>- THE BATTERIES ARE NOT INSTALLED CORRECTLY</li> <li>- BATTERIES ARE DEAD</li> </ul>	<ul style="list-style-type: none"> <li>- FIT BATTERIES</li> <li>- CHECK THE POLARITY</li> <li>- REPLACE THE BATTERIES.</li> </ul>
THE SCREEN SWITCHES ON BUT STIMULATION IS AT 0	<ul style="list-style-type: none"> <li>- THE ELECTRODES ARE NOT IN CONTACT WITH THE SKIN</li> <li>- INTENSITY LEVEL AT 0</li> <li>- SOFTWARE FROZEN</li> <li>- LOW BATTERY POWER</li> </ul>	<ul style="list-style-type: none"> <li>- APPLY GEL AND PLACE ELECTRODES IN CONTACT WITH THE SKIN</li> <li>- INCREASE THE INTENSITY</li> <li>- RESET THE APPLIANCE BY REMOVING THE BATTERIES REFIT THE BATTERIES AFTER 1 MINUTE</li> <li>- REPLACE THE BATTERIES.</li> </ul>
THE SCREEN SWITCHES ON BUT STIMULATION IS NOT NOTICEABLE	<ul style="list-style-type: none"> <li>- LOW INTENSITY</li> <li>- NOT ENOUGH GEL ON THE ELECTRODES</li> <li>- GEL DIFFERENT TO THE ONE THAT COMES WITH THE PRODUCT USED</li> <li>- THE BELT HAS NOT BEEN FITTED CORRECTLY</li> <li>- THE ELECTRODES ARE DIRTY</li> <li>- LOW BATTERY POWER</li> </ul>	<ul style="list-style-type: none"> <li>- INCREASE THE INTENSITY</li> <li>- APPLY MORE GEL</li> <li>- APPLY THE GEL THAT COMES WITH THE PRODUCT</li> <li>- FIT THE BELT CORRECTLY</li> <li>- CLEAN THE ELECTRODES</li> <li>- REPLACE THE BATTERIES.</li> </ul>
THE BUTTONS DO NOT FUNCTION CORRECTLY	<ul style="list-style-type: none"> <li>- SOFTWARE FROZEN</li> <li>- LOW BATTERY POWER</li> </ul>	<ul style="list-style-type: none"> <li>- RESET THE APPLIANCE BY REMOVING THE BATTERIES REFIT THE BATTERIES AFTER 1 MINUTE</li> <li>- REPLACE THE BATTERIES.</li> </ul>

## MAINTENANCE

1. Clean the electrodes after using a soft cloth lightly moistened with water.
2. Make sure that no water gets inside the appliance.
3. Remove the batteries from the unit when it will not be used for prolonged periods. Leaking batteries may damage the unit.
4. Used batteries must be deposited in special collection containers for toxic waste. The law requires proper disposal of batteries.
5. Store the appliance in a dry, well-ventilated place.

## TECHNICAL SPECIFICATIONS

MODEL	e-SHORT YR40
FREQUENCY	1Hz – 60Hz
BATTERIES	2 X 1.5V AAA ALKALINE BATTERIES
ELECTRODE SURFACE	546 Cm <sup>2</sup>
E-SHORT WEIGHT	400 Gr

## Warranty

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25<sup>th</sup> May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10<sup>th</sup> July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT n<sup>o</sup> B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

<b>Date of purchase.-</b>	<b>Purchased in.-</b>
<b>Purchaser.-</b>	<b>Dealer stamp.-</b>

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.