# YU110 SONIXVISAGE

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[www.tecnovita.eu](http://www.tecnovita.eu)

# Welcome

We would like to thank you for purchasing our product and for the trust you have placed in Tecnovita.

All of our products have been tested and submitted to strict quality controls.

Read this manual before use the product.

Get professional results anywhere due to its small size and high quality

Due to the advance technology developed in the SonixVisage you can show off a mositurise, firmly and elastic skin.

Its small size and recargable battery achieves you to enjoy an ultrasound threatment anywhere you need.

## Benefits

Softens wrinkles and expression lines.

Corretcs crow´s feets

Eliminates puffiness and dark circlces

Bar code lines

Combat saggy skin

Reduce jowls

Stimulate the production of caollagen

Eliminate dead skin cells.

Reduce the skin fat to get an esplendid cell.

# Ultrasonic Use

## **How ultrasonic works?**

The high frequency ultrasonic wave technology helps to improve the cell membrane metabolism, contributing to better blood circulation. This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat both the sagging skin and also remove wrinkles, dark circles, bags under the eyes, etc..

The ionic technology improves facial cleansing, helping to remove grease, dirt and traces of makeup.

## **What are the advantages over other methods?**

3MHz ultrasonic waves are safe, painless and non-invasive and they also aid sonophoresis, allowing the nutrients in skin care creams to penetrate deeper.

## **Is it painful?**

It is completely painless and safe.

# Ultrasonic Effects

## Mecanic Effect


Ultrasonic technology offers you a cellular micro massage to increase skin permeability. This allows to increase the permeability of membranes, which helps for a penetration at the skin

## Thermal Effect

The mechanical energy absorbed by the skin is transformed into thermal energy, which in turn generates heat. This process helps the stimulation of cell metabolism and circulation.

## Chemical Effect

The chemical action of the ultrasound activates oxygen molecules, creating better oxygenation of the skin and therefore better penetration therein of cosmetic products.
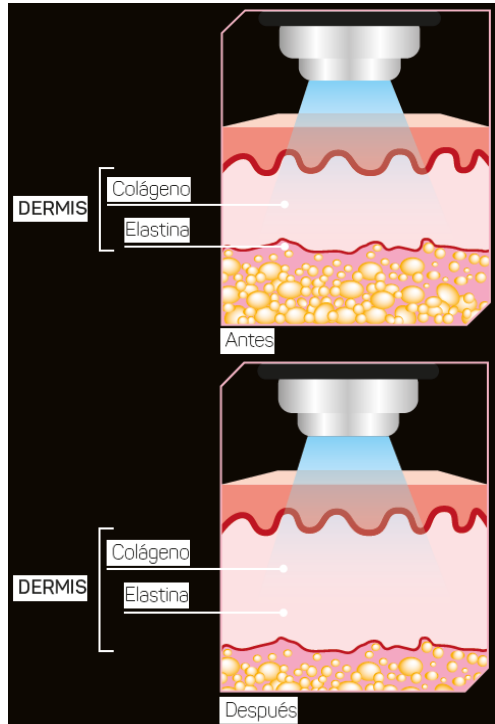
# LED USE

## How Led technology works in our skin?

Led technology is practically new. It helps us to rejuvenate our skin with the use of different colors and different wavelengths. Depending on the results that we want, we have to choose between different colors.

## Red LED

The red LED light stimulates the production of collagen, helping the production of collagen and elastin, obtaining a firm and SUPPLE SKIN. Reduce la apariencia de arrugas y ayuda a retrasar su desarrollo. The red LED light stimulates the blood circulation. This process makes new cells and makes the skin look healthier and rejuvenated.

## Blue LED

Blue light have purifying and antibacterians properties. It's specially indicated for fat and sensible skins.

## Green LED

It's useful for pigmented skin, skin blemish and sunspots. It also have calming and relaxing properties.

## Technical characteristics

Advance technology with a stylish design.

Ultrasonic facial function.

3 Intesity levels in ultrasonic function.

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3 LED colors.

Rechargeable massage head.

Great confort due to their small size.

# Specifications

**Output power:** 3 MHz

**Weight:** 190 gr.

**Mains adapter.** 100-240 VAC 50-60 Hz 0.2Amax

**Dimensions:** L x An x Al (4'5 x 4'7 x 16'5 cm).

## SonixVisage operating instructions

**Precautions: Consult your doctor before beginning to use this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.**

- This appliance has been designed for home use.
- This appliance has been designed only for use on the face.
- **Do not use the head on your throat, neck, eyeballs, ears or ear canals.**
- We recommend working on clean skin, without makeup, so that the ultrasound penetrates better.
- Before use remove all accessories such as necklaces, earrings, contact lenses, body piercings, etc..
- In the ultrasound mode never use the appliance without first applying plenty of ultrasound gel. If you do, then the ultrasound will not be transmitted correctly and it can also damage the appliance.
- In the ultrasound mode do not stop moving the appliance. Never let it rest on the skin as this could cause burns.
- In the ultrasound mode do not use the appliance for more than 20 minutes.
- In the ultrasound mode use always the protection glasses.
- In the LED mode do not use the appliance for more than 15 minutes.
- Do not use the appliance if the cable is damaged.
- Do not unplug the appliance by tugging on the cable.
- Do not move the appliance by pulling on the mains cable.
- Do not touch the appliance with wet hands
- Always clean the head after use, using a cloth slightly dampened with water. Never immerse the appliance in water.
- Do not expose the unit to high temperatures or magnetic fields.
- If you notice that the appliance is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

# Package contents

1. Ultrasound gel YUG01
2. YU110 SONIXVISAGE
3. Charge base
4. Mains adapter
5. Charge cable
6. Protection glasses

# Quick Start


## Charging the battery

- Fully charge the battery before using the appliance for the first time. You can charge it through the base or directly from the power adapter
- It will take approximately 3 hours to fully charge the battery, although the first charge could take up to 16 hours.
- The LED (B) will light up red when the battery is low.
- When the battery is completely dead the appliance will beep three times and switch off automatically.

To prolong battery life always recharge the battery after it is fully discharged and do not leave the battery on charge longer than necessary to charge it completely.

## Charge base

1.

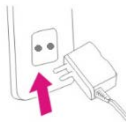
2.

## Charge cable

1.

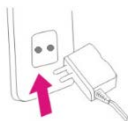2.

## Control Functions

- |                 |                                    |
|-----------------|------------------------------------|
| 1. ON/OFF       | Encendido/Apagado                  |
| 2. LED          | Active/Select LED mode             |
| 3. ULTRASOUND   | Active/Select Ultrasound intensity |
| 4. CHARGE       | Charge indicator                   |
| 5. POWER        | Power indicator                    |
| 6. Level LOW    | Low intensity ultrasound           |
| 7. Level MIDDLE | Middle intensity ultrasound        |
| 8. Level HIGH   | High intensity ultrasound          |

## Beginning a session

- Fully charge the battery before using the appliance for the first time.
- Before beginning a session, remember to remove makeup and all accessories such as necklaces, earrings, contact lenses, piercings, etc... Make sure your skin is clean and dry.
- In ultrasound mode always apply an abundant amount of gel to the area and begin to work it by slowly sliding the head over the area but without pressing, as shown in the diagrams contained in the Application Zones section of this manual.

Before starting the LED mode remember to put on the protection glasses

If you need to purchase ultrasound gel, you can find it in the shop where you purchased SonixVisage or at our website: [www.tecnovita.eu](http://www.tecnovita.eu)

# For correct, safe use of SonixVisage

For correct and safe use of SonixVisage, you should keep the following application times in mind:

- **Per work area:** By "work area" we mean each part of the face where SonixVisage can be applied. The various work areas are:

- o *Forehead:* **max. 3 minutes**
- o *Around the eyes:* **max. 3 minutes**
- o *Cheeks and cheekbones:* **max. 3 minutes**
- o *Around the lips:* **max. 3 minutes**
- o *Chin and jowls:* **max. 3 minutes**

- **Per daily session: 15 minutes maximum.**

By session we mean the application of SonixVisage on different work areas in the same day.

- **Between sessions: 48 hours minimum.**

By "between sessions" we mean the rest time between each session.

*E.g.: In one week, a session schedule would be: Monday, Wednesday, Friday and Sunday.*

- **Every 24 sessions (approx. 1 and a half months of treatment): rest 1 and a half months.**

After this rest period you can repeat the treatment.

## **VERY IMPORTANT**

**\*\*\* Never exceed the application times (a total of 15 minutes per session in one day). Exceeding recommended applications can cause fatigue and skin irritation.** Only short applications should be applied to allergic or sensitive skin.

**\*\*\* And, above all, respect the rest period of 48 hours between sessions. This is the time that your skin needs to rest.**


# Areas of application

## Forehead

*Apply gel over the area to be worked*

Slide the handset slowly, without pressing, in an upward circular motion, covering the entire area of the forehead, as shown in the picture.

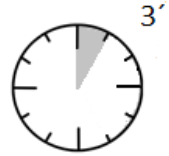
## Around the eyes

*Apply gel over the area to be worked*

The massage movement should be applied from inside out away from the eye area as shown in the picture. It is important to massage in the correct direction as this will produce better results.

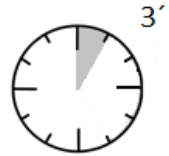
\* Remember **DO NOT** apply *SonixVisage* on the eyeball or upper eyelid.

## Cheeks and cheekbones

*Apply gel over the area to be worked*

Slide the handset, using circular movements, from the start to the end of the cheekbones and cheeks covering the entire area.

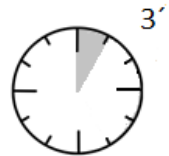
\* Remember **DO NOT** apply *SonixVisage* on the ears.

## Around the lips

*Apply gel over the area to be worked*

Perform gentle movements from the centre outward around the whole lip area, as shown in the picture.

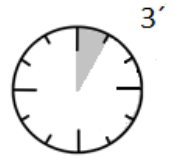
## Chin and jowls

*Apply gel over the area to be worked*

Slide the handset using an upward circular motion from the centre outward, covering the entire area.

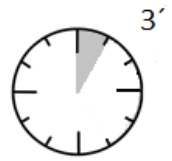
***\*Remember DO NOT apply SonixVisage on the neck or throat***

## LED Mode

- Before starting the session LED remember to put on protection glasses.
- Do not apply over the same point 5".
- Do not apply more than 10 'for color.

## Luz Roja

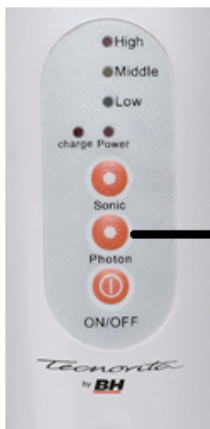x1

Rojo  
Red  
Rouge  
Vermelho  
Rosso  
Rot

## Luz Azul

x2

Azul  
Blue  
Bleu  
Azul  
Blu  
Blau

# Luz Verde

x3

Verde  
Green  
Vert  
Verde  
Grün

## Frequently asked questions about SonixVisage

### What is SonixVisage?

The high frequency ultrasonic wave technology helps to improve the cell membrane metabolism, contributing to better blood circulation. This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat both the sagging skin and also remove wrinkles, dark circles, bags under the eyes, etc..

### What are the advantages over other methods?

3MHz ultrasonic waves are safe, painless and non-invasive and they also aid sonophoresis, allowing the nutrients in skin care creams to penetrate deeper.

### Is it painful?

It is completely painless and safe.

### Is it normal to hear a slight buzzing in your ears while using SonixVisage?

Yes, some users may detect a slight buzzing sound during its application.

### Can it be used on all parts of the body?

The 3Mhz SonixVisage technology has been specifically designed for use on the face. **Under no circumstances must SonixVisage be applied to the throat, neck, eyeballs, ears and ear canals.**

**Is SonixVisage effective on “Crow’s feet”?**

Yes, SonixVisage is effective on "crow's feet", aiding the regeneration and tighten of wrinkle-causing cells.

**Can SonixVisage be used to remove wrinkles?**

Yes, SonixVisage is an effective way to remove wrinkles.

**After how many sessions will I begin to see results?**

You should notice results after the first five sessions, however results may vary depending on the skin of each individual.

**How long should I wait between sessions?**

It is very important to respect the 48 hour break between sessions, as this is the time that your skin needs to rest.

**Do wrinkles reappear after using SonixVisage?**

Although SonixVisage ensures effective and long-lasting results, wrinkles are part of the natural ageing process. To reduce and delay the onset of wrinkles SonixVisage maintenance sessions are recommended.

**Are there any side effects?**

If you are healthy, SonixVisage has no reported side effects, nevertheless in some people it may produce a reddening of the skin that can last for a few minutes.

**What it ultrasound gel for? When is it used?**

The ultrasound gel is essential, ALWAYS use is before beginning a session. This will ensure that the ultrasonic waves penetrate correctly and prevent damage to the unit.

**Can I use SonixVisage without applying gel?**

Never, because the gel acts as a conductor for the ultrasonic waves on your face. If you do not use gel then the SonixVisage application loses its effectiveness and the machine can be severely damaged.

**Can I use any type of gel?**

We recommend using only TECNOVITA gel, as the use of any other type may damage the machine.

**Is it necessary to drink water?**

No, it is not strictly necessary but it is advisable to hydrate yourself. The recommended quantity is at least 2 litres during treatment, as this helps to hydrate the skin.

### **When is it inadvisable to use SonixVisage?**

- if you suffer from hyperpigmentation, broken capillaries
- if you have sunburn
- if you suffer from severe skin problems
- active skin infections
- facial paralysis
- If you are pregnant or think you might be
- If you are breast feeding.
- undergoing fertility treatments.
- couperosis
- temporal arthritis
- malignant tumours
- diabetes
- kidney problems
- physical or mental deficiencies
- metal implants or prostheses in the area of application
- scarring in the area of application
- physical malformations.
- epilepsy
- long-term treatment with medications: anti-inflammatories, anticoagulants, antibiotics or antihistamines.
- treatment for psychological or psychiatric disorders

### **Can I use SonixVisage if I am pregnant?**

NEVER use SonixVisage if you are pregnant or suspect that you might be.

### **Can SonixVisage be used during lactation?**

Do not use it if breastfeeding.

### **Can SonixVisage be used on the breasts?**

Never apply it to the breasts.

### **Can SonixVisage be used with skin conditions?**

No, It should not be used if you have any serious skin condition. Consult you doctor.

### **Can I use SonixVisage if I have hearing implants?**

No, it should not be used if you have hearing implants fitted.

### **Can I use SonixVisage if I have implants in my mouth, gold teeth, metal implants or similar?**

Do not use it if you have these types of implant, as it may cause discomfort or tenderness in the area of application.

### **Can I use SonixVisage when I am bathing or showering?**

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

### **How do I clean SonixVisage?**

Never immerse SonixVisage in water to clean it. Use a cloth slightly dampened with water.

## **Warranty**

We suggest that you read through this document carefully in order to fully understand the cover that the manufacturer offers as well as your responsibilities as owner.

### **GUARANTEES**

#### **GENERAL INFORMATION**

EXERCYCLE, S.L. would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read through these warranty conditions and then get in touch with your local authorised dealer.

This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25<sup>th</sup> May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10<sup>th</sup> July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

#### **1. GUARANTEE PERIOD**

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later.

It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries.

The stated warranty period is only applicable to Europe.

## 2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.

- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.

- Repairs will be free for the consumer with regard to labour, materials and carriage.

- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended temporarily until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2 year period be initiated.

The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer. The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the

contract (ii) will not proceed when the non-conformity is of minor importance.

### 3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).
- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from lack of maintenance (accumulation of dirt on any part of the product or its accessories, etc.). (The User Manual provides advice on how to install and use the product).
- Repairs carried out by persons not authorised by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorised by the manufacturer.
- Deviations from the standard product, such as light vibrations, noises or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks, falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

### 4. CLAIMING AGAINST THE GUARANTEE

- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.
- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
- This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
- With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law.

Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

## **Certificates**

By virtue of this document, Exercycle S.L. declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC

Please keep the receipt as proof of purchase.

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