*Tecnovita*  
by **BH**

**Instrucciones de uso y garantía**  
**Operating instructions and guarantee**  
**Mode d'emploi et garantie**  
**Instruções de utilização e garantia**  
**Manuale d'uso e garanzia**  
**Bedienungsanleitung und Garantie**

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**902 170 258**

Thank you for purchasing our YM301 Feet Care Presotherapy. All of our products undergo rigorous testing to ensure maximum quality. To use the Massager correctly, first read through the whole manual before using it.

## INTRODUCTION

The YM301 Feet Care Presotherapy massager uses the traditional Shiatsu and acupressure massage technique. Thanks to the latest state-of-the-art technology you can enjoy a real massage that helps to relax your muscles, eliminate fatigue, combat stress and tension. The device also features a heat function.

## GENERAL SAFETY INSTRUCTIONS.

1. Inspect the device before use to ensure that it is in good condition and that the cable is not damaged. If you notice a defect or missing parts, contact the distributor of the product.
2. If you have any health problem, consult your medical doctor before using the device.
3. The device must not be used by children or the disabled without suitable adult supervision.
4. Excessive use can make the massage device overheating. If this happens let the device get cool.
5. If you experience any discomfort while using the unit, cease use immediately and consult your doctor.
6. To avoid burns, use the heat function with caution. Do not use it if you have a skin problem (sensitive skin, poor circulation, ...).
7. Do not cover the unit during operation (e.g. with a blanket or pillow) because it may overheat and cause damage.
8. Do not use the unit while sleeping.
9. This device is designed to be used on the feet. Do not use it on any other part of the body.
10. Do not use the unit for longer than 30 minutes at a time.
11. Do not use the unit just before or after eating. Allow at least 1 hour.
12. Do not use the device in damp/humid surroundings and prevent it from getting wet. Do not touch the device with wet hands while it is operating.
13. Do not use the unit if the ambient temperature is above 40°C.
14. Use the unit under normal room temperature conditions. If you take the unit out of


surroundings at a low temperature allow at least one hour before turning the device on as internal condensation might occur and affect how the unit operates.

15. Do not tighten or force any body part in the massage mechanism in motion.
16. Do not insert objects inside the unit.
17. Do not stand on the device.
18. Do not grip the appliance by the cable, do not crush it
19. Do not move the device by pulling the cable.
20. To avoid breakage, do not wrap the cord around the device.
21. Switch off the device completely if the appliance is not going to be used for a while and unplug it.
22. Do not unplug the appliance without switching it off first.
23. In the event of a power cut disconnect the plug from the mains supply.
24. To disconnect the device from the main power supply, pull the plug and not the cable.
25. To prevent the risk of fire, damage or electric shock, use the correct voltage
26. Do not expose the device to high temperatures. Keep it in a dry place.
27. If you notice that the device is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty. Contact your technical assistance service.

Consult your doctor about using the device in the following cases:

1. You are on medical treatment
2. Hypertension, heart disease, arteriosclerosis, cerebral thrombosis
3. Malignant tumours
4. Pregnancy
5. Osteoporosis
6. Using a pacemaker or other internal device
7. Acute diseases

## DESCRIPTION

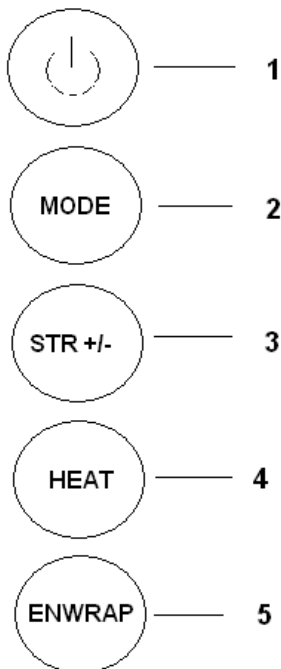- 1 - On/Off button
- 2 - Mode selection - 2 available modes, red or blue light
- 3 - Massage intensity setting - 3 levels - Red light for high, blue light for low level
- 4 - Select Heat Mode - Lights up red if the mode is enabled (enabled by default)
- 5 - Air massage mode

## **USE**

1. Place the device on a flat surface.
2. Connect the device to the main power supply.
3. Press the on/off (1) button to switch the appliance on. The default settings for the device are heat function on (red light), high intensity level (red light), mode 1 (light blue).
4. To switch the heat function off press the HEAT (4) button, the red light will switch off.
5. To change the intensity level, press the STR (3) button. Blue light – low level, red light – high level.
6. To change the mode, press the MODE (2) button. Blue light – mode 1, red light – mode 2.

7. To receive only air pressure massage, press the ENWRAP (5) button.
8. To switch the device off, press the on/off (1) button.

## **MAINTENANCE AND CLEANING.**

1. Unplug the device and leave it cool before cleaning it.
2. Clean it with a soft dry cloth. Do not allow the appliance to come into contact with water or any other liquid. Do not immerse the appliance in liquid to clean it. Do not use abrasive cleaners, brushes, gasoline, glass or furniture polishes or solvents to clean it.

## TECHNICAL SPECIFICATIONS

MODEL	YM301 FEET CARE
VOLTAGE	DC24V/2A
POTENCIA	45

## Warranty

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25<sup>th</sup> May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10<sup>th</sup> July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

<b>Date of purchase.-</b>	<b>Purchased in.-</b>
<b>Purchaser.-</b>	<b>Dealer stamp.-</b>

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.