*Tecnovita*

by **BH**

**Instrucciones de uso y garantía**  
**Instructions for use and warranty**  
**Instructions pour l'utilisation et la garantie**  
**Instruções de uso e garantia**  
**Bedienungsanleitung und Garantie**  
**Manuale d'uso e garanzia**

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**902 170 258**

Thank you for purchasing our AURA YM851 massage helmet.

All of our products undergo rigorous testing to ensure maximum quality.

To use the AURA YM851 massage helmet correctly, first read through the manual before using it.

### **INTRODUCTION.-**

The AURA YM851 massage helmet is a unit that combines massage around the temple area with massage around the occipital and neck area.

The AURA YM851 massage helmet operates using a microchip control system, applying the modern technologies of air pressure and vibration which, along with other features such as heat and music, provide a sensation of total massage, eliminating stress and tension and stimulating brain activity.

Specific uses:

- Anti-stress therapy eliminates stress and fatigue
- Combats facial tiredness (rings and bags under the eyes, wrinkles).
- Provides a feeling of wellbeing and mental massage.

### **DESCRIPTION.-**

- The product is a multi-functional high-tech device for massaging and relaxing the eye, temple, occipital and neck area.
- Thanks to an advanced microchip, this product combines different variables (vibration, air pressure, infrared heat and relaxing sounds) with full control via a remote control unit.
- The helmet AURA YM851 produces an enveloping vibratory motion that stimulates blood flow by means of various multi-frequencies.
- The helmet AURA YM851 uses an infrared based heat emitter that activates the radiated cells, causing the molecules to vibrate, strengthening their bonding properties and improving the regenerative capacity of the skin tissue.
- The product successfully applies the concept of air pressure.

### **SAFETY INSTRUCTIONS.-**

1. Inspect the device before use to ensure that it is in good condition and that the cables are not damaged. If you notice a defect or missing parts, contact the distributor of the product.
2. Do not use this product together with other therapeutics devices.
3. The device must not be used by children or the disabled without suitable adult supervision.
4. This device must not be used as a replacement of any kind of treatment.
5. Do not use this device under the effect of alcohol and/or drugs.
6. Please avoid falling asleep while using this device
7. Please do not use this massager immediately after meal.
8. Do not drive while using the device.
9. Do not lie down while using this device
10. Please adjust the sound volume at a reasonable level, high volume may lead to hearing problems.
11. To prevent the air chamber from over-pressurising, first place the helmet on and then start the device.
12. Do not use the massage device in damp/humid surroundings and prevent it from getting wet. Do not touch the device with wet hands while it is operating. Under no circumstances must this product be used in the bath or shower.
13. We recommend to keep your eyes closed for a few seconds and then opening them slowly after a massage session.
14. Do not expose the device to direct sunlight or corrosive gases.
15. Switch the device off when you have finished using it.
16. Switch the device off before cleaning it.

**17.** Do not wash the product under a tap or immerse it in water. Use a damp cloth, making sure to dry the device thoroughly before using it again.

**18.** Do not expose the device to high temperatures. Store it in a dry place.

**19.** If you notice that the AURA YM851 massage helmet is not operating correctly, do not attempt to repair, dismantle or modify it. Unauthorised repairs may cause injury, in addition to breaking the terms of the warranty.

Contact your technical assistance service or local distributor.

**20.** This device has been designed for domestic use only not commercial.

## CONTENT OF THE PACKAGE


1. MASSAGE HELMET
2. REMOTE CONTROL
3. CHARGER ADAPTOR
4. USB CABLE
5. EARPHONES
6. TF MEMERY CARD

## DESCRIPTION

1. Height regulation knob
2. Massage nails
3. Eye massager
4. Widht ajustement knob
5. Master unit
6. Switch button
7. Features screen
8. Programs selections button
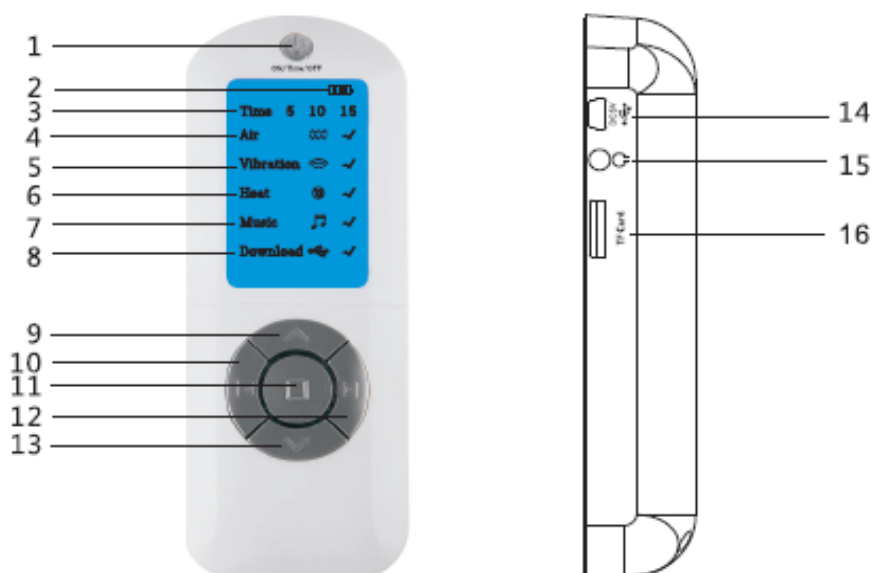
## REMOTE CONTROL

1. Power button (Long press of 2 sec to turn on/off the power. Short press to set the time).
2. Battery level indicator.
3. Time indicator – 5, 10 and 15 min.
4. Aire pressure indicator.
5. Vibration indicator.
6. Heating indicator.
7. Mp3 player indicator
8. Mp3 download indicator
9. Upward selection button
10. Modes or music adjustment button
11. Massage mode (automatic or free) selection button and music stop
12. Modes or music adjustment button
13. Downward selection button
14. Power charger and Mp3 download interface
15. Earphones interface
16. Memory card slot


## BEFORE USE

Before starting using the massage helmet YM851 AURA, please make sure that the device is fully charged. Wear the helmet and adjust its size to the size of the head through the knobs 1 and 4.

After adjusting the helmet, press the “” button to switch on the device. The helmet will start operating automatically in a preset program. To switch off the device press the “” button again during few seconds.


## AUTOMATIC MASSAGE MODE

Under Auto Massage mode, press  button to choose between the 6 following different automatic Massage programs.

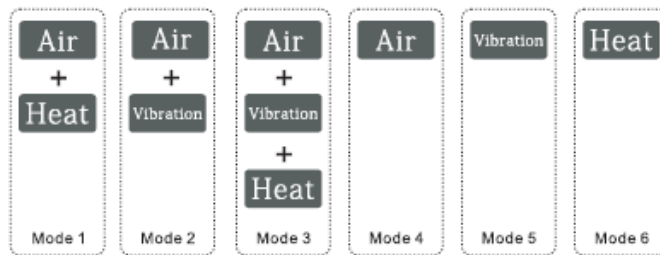
## CUSTOM MESSAGE

### AIR

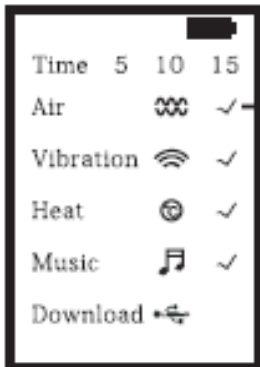Move up or down (button 9 and 13) to the air pressure selection bar. The selection will blink. You can use the left or right button (10 and 12) to adjust the 4 submodels of air pressure.

### VIBRATION

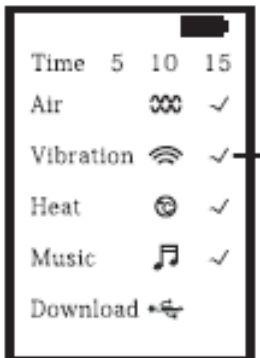Move up or down (button 9 and 13) to the vibration selection bar. The selection will blink. You can use the left or right button (10 and 12) to adjust the 10 submodels of vibration.

### HEAT

Move up or down (button 9 and 13) to the heating selection bar. The selection will blink. You can use the left or right button (10 and 12) to adjust the two submodels of heating (38°C and 42°C).

## MUSIC FUNCTION

- This device is equipped with the Autoplay function. When the device is turned on, music will start automatically. On the remote control, the music option will be marked. To stop music, press the  button. To turn on music again, press the  button again. To choose the song to be played use buttons 10 and 12.
- To download music into the massage helmet, connect it to a pc through the output n° 14 and switch it on with the “” button. The pc will detect the TF memory card automatically.

**NOTICE:** this device only supports Mp3 Format files.

## MESSAGE TIME ADJUSTEMENT

This device has 3 preset massage durations: 5, 10 and 15 minutes. To change it, press the “” until the function time blinks on the remote control.

## MAINTANANCE, CLEANING AND STORAGE

1. Switch off the appliance and allow it to cool down before cleaning it.
2. Clean it with a slightly dampened sponge or soft cloth or with cotton and alcohol.
3. Do not allow the appliance to come into contact with water or any other liquid. Do not immerse the appliance in liquid to clean it. Do not use abrasive cleaners, brushes, gasoline, glass or furniture polishes or solvents to clean it.
4. Keep the device under dry and ventilated conditons, away from direct sunlight and corrosive gases.

## HOW TO DISPOSE ELECTRICAL DEVICES.-

1. If a product bears the crossed out wheeled bin symbol, this means that European Directive 2002/96/EC applies.
2. Waste electrical or electronic equipment must not be disposed of in household waste, but at the collection points designated by the town Council or local authorities.
3. Proper collection and processing of unusable devices helps avoid potential risks to the environment and public health.
4. For more information about how to dispose of waste electrical and electronic equipment, contact your town hall, waste disposal service or the establishment where you purchased the product.

**Keep these instructions safe.**

## TECHNICAL SPECIFICATIONS

MODEL	AURA YM851
ENVIRONMENT OPERATING TEMPERATURE	5°C – 40°C
RELATIVE HUMIDITY	< 80%
ELECTRICAL SUPPLY	5V-DC
POWER	5W
CONTROL UNIT DIMENSIONS	15 x 6 x 3 cm.
DIMENSIONS, MAIN UNIT	23x27,8x24,3cm.
NET WEIGHT	1974g.
POWER ADAPTOR	AC100-240V DC5V 1000MA

**THE MANUFACTURER RESERVES THE RIGHT TO MODIFY THE SPECIFICATIONS OF ITS PRODUCTS WITHOUT PRIOR NOTICE.**

### Garantía

This product is guaranteed for 2 years in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25<sup>th</sup> May, 1999, regarding certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10<sup>th</sup> July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product.

<b>Date of purchase.-</b>	<b>Purchased in.-</b>
<b>Purchaser.-</b>	<b>Dealer stamp.-</b>

Please keep the receipt as proof of purchase.

Should you have any queries please do not hesitate to phone our customer service line:

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BH China Co., Ltd.