ESPAÑOL/ENGLISH

MANUAL DE USUARIO
USER MANUAL

ESPAÑOL

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ENGLISH

We would like to thank you for purchasing **SonixPlus** and for the trust you have placed in **Tecnovita by BH**. By means of this handy User Manual and following its simple instructions you will achieve results comparable to those obtained from the hands of a professional beautician.

We strongly recommend that you read this User Manual thoroughly. Reading and understanding it will help you use **SonixPlus** effectively and safely.

Contents

Chapter 1

- Contents of the box

Chapter 2

- Basic characteristics of SonixPlus
- SonixPlus operating instructions

Chapter 3

- SonixPlus application modes

Chapter 4

- For correct, safe use of SonixPlus

Chapter 5

- Frequently asked questions

Chapter 6

- Important safety instructions
- Warranty
- Certificates
- Technical specifications
- Parts list

Chapter 1

Contents of the box

Inside this box you will find everything you need to enjoy SonixPlus from now onward. Carefully check the following illustration showing all of the elements that the pack should contain.

1 – Main unit

2 - Flat handset

3 - Tapered handset

4- Ultrasound gel

5 – User manual


Chapter 2

Basic characteristics of SonixPlus

SonixPlus is a hi-tech device designed to stimulate skin tissue and improve the elasticity of the skin.

The technology used employs the emission of high frequency ultrasonic waves which help to improve the cell membrane metabolism, contributing to better blood circulation. This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat both the sagging skin and also remove wrinkles, dark circles, bags under the eyes, etc.

Characteristics for model SonixPlus YU180

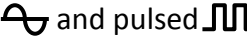- **Time:** Programmable from 1 to 60 minutes.
- **Energy:** Programmable at 10 intensity levels
- **Waveform mode:** Selection of 2 wave types: continuous and pulsed 
- **Handle mode:** Selectable: left (L) or right (R)

There are only three simple functions needed to program this unit: time, energy and wave type, in addition to selecting the handset that you wish to use according to the area of the face to be treated.

The **time** function allows you to choose how long you wish to work with SonixPlus. You can reduce the maximum time of 60 minutes by 1 minute intervals. Once the selected time has elapsed, SonixPlus will inform you with a beep.

The **energy** function controls how deep the sound waves penetrate your skin. The greater the intensity, the deeper the waves reach.

SonixPlus also lets you choose the wave type, either **continuous** or **pulsed**

helping you select the time interval during which the waves reach your skin.

In addition, you can independently select either of the **two handsets**, one with a flat base and the other with a tapered end, depending on the area where they will be applied. The different sizes and shapes ensure that they fit perfectly into tricky areas of the face such as around the contours of the eyes and lips. You must select one or the other, it is not possible to work with the 2 handsets at the same time.

SonixPlus operating instructions

Precautions: Consult your doctor before beginning to use this appliance. This warning is especially important for people with health problems. Use by those under 18 years of age is not recommended. Keep this appliance out of children's reach.

- We recommend working on clean skin, without makeup, so that the ultrasound penetrates better.
- We recommend removing all accessories such as necklaces, earrings, contact lenses, piercings, etc..
- Never use the appliance without first applying plenty of ultrasound gel. If you do, then the ultrasound will not be transmitted correctly and it can also damage the appliance.
- Do not stop moving the handset. Never let it rest on the skin as this could cause burns.
- **Do not use the handset on your throat, neck, eyeballs, ears and ear canals.**
- Always clean the handsets using a cloth slightly dampened with water. Never immerse them in water.

Starting

- Connect the **handsets** to the **main unit**. To do this, take the plugs for both **handsets** and insert them into the sockets at the back part of the **main unit**. Each of the handsets can be plugged into either of the two connections.
- Take the **mains cable** and plug it into the electric wall socket (220V).

Monitor elements

The SonixPlus monitor consists of seven keys, which have the following functions:

- **TIME + - (6, 7)**

Adjustable from 1 to 60 minutes (in steps of 1 minute)

- **ENERGY + - (4, 5)**

Divided into 10 levels: 1-2-3-4 low intensity (green light), 5-6-7 medium intensity (yellow light) and 8-9-10 high intensity (red light).

- **WAVEFORM MODE (2)**

Adjustable for continuous  or pulsed  emission.


- **HANDLE MODE (3)**

To select the desired handset, depending on the type of application (flat or tapered).

Left (L) or right (R).

- **START (1)**

The red light indicates that it is switched on.

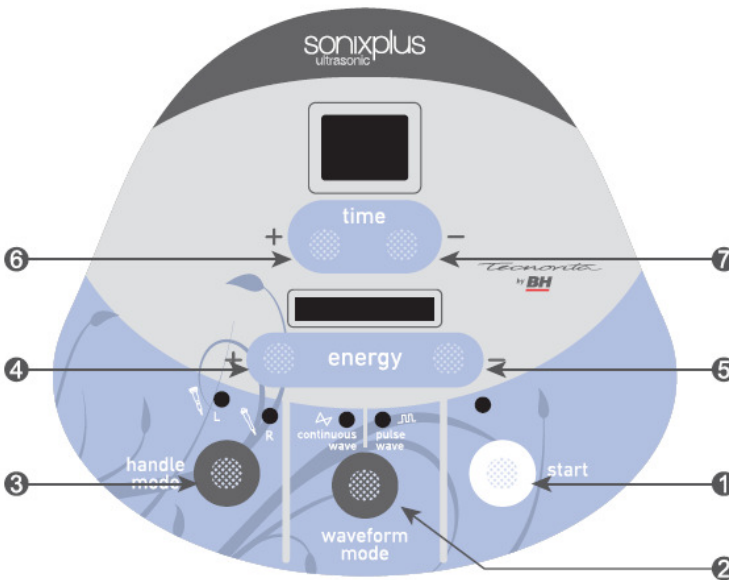Beginning a session

Once you have plugged the machine into the mains supply press the start (1) button to switch the unit on.

The following settings are established by default:

TIME: 30

WAVEFORM MODE: continuous 

HANDLE MODE: left (L)

ENERGY: 0

From here you can change the time, energy and wave type settings in order to establish the type of session you want:

Changing the duration of the session (*TIME*)

You can change the time setting using the + and - keys (6), (7)

The displayed default value is 30 minutes but this can be changed in intervals of 1 minute between values of 1 and 60.

Changing massage intensity (*ENERGY*)

To select the energy level, press the + and - keys (4), (5).

You can select from 4 low levels (green light), 3 medium (yellow light) or 3 high (red light).

Changing work mode (*WAVEFORM MODE*)

You can change to whichever work mode you want for your session by pressing **waveform mode** (2) choosing between continuous or pulsed mode.

Selecting the desired handset (*HANDLE MODE*)

You can select whichever handset you want for your session by pressing **handle mode** (3). Choose between the left handset (L) or the right handset (R).

Finally, take hold of the selected handset to begin your session. Apply the ultrasound gel to the area of the face that you wish to work and move the handset over it as indicated in the “Application modes” section.

You will see the time countdown in steps of 1 minute

You can pause stop the session by pressing the button start (1), the unit will turn off. By pressing the start button (1) you will return directly to the previous session.

When the time counter reaches **00** the unit will stop and emit a beep to tell you that the selected session has finished and it will turn off.

Chapter 3

SonixPlus application modes


Before beginning a session, remember to remove all accessories such as necklaces, earrings, contact lenses, piercings, etc.

Always apply an abundant amount of gel to the area and begin to work it by slowly sliding over it without pressing, as shown in the diagrams.

If you need to purchase gel or any of our preparation creams you can find them in the shop where you purchased SonixPlus or at our website: www.tecnovita.eu

Areas of application

Forehead

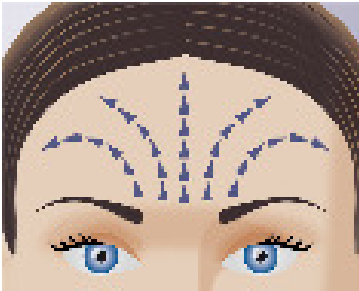
Apply gel over the area to be worked

Slide the handset slowly, without pressing, in an upward circular motion, covering the entire area of the forehead, as shown in the picture.

Max. Time	Energy	Type
3 min.	Medium-high (5 to 10)	continuous 

Eye contour

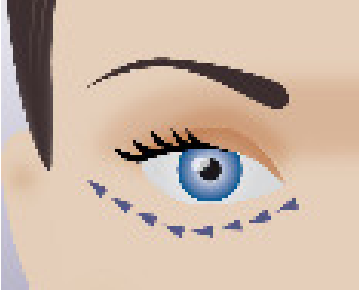
We recommend using the tapered handset for a more accurate application in the area around the eyes.

**** Remember DO NOT apply SonixPlus on the eyeball or upper eyelid.***

The massage movement should be applied from inside out away from the eye area as shown in the picture. It is important to massage in the correct direction as this will produce better results.

Max. Time	Energy	Type
3 min.	Low (1 to 4)	pulsed 

Cheeks and cheekbones

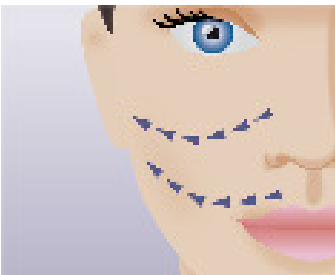
Apply gel over the area to be worked.

Slide the handset, using circular movements, from the start to the end of the cheekbones and cheeks covering the entire area.

**** Remember DO NOT apply SonixPlus on the ears.***

Max. Time	Energy	Type
3 min.	low-medium (1 to 7)	pulsed 

Lip contour

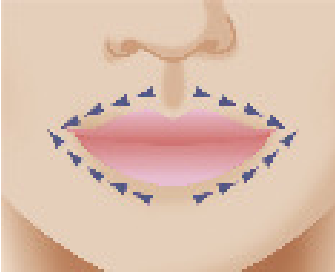
We recommend using the tapered handset for a more accurate application in the area around the lips.

Perform gentle movements from the centre outward around the whole lip area, as shown in the picture.

Max. Time	Energy	Type
3 min.	Low (1 to 4)	pulsed 

Chin and jowls

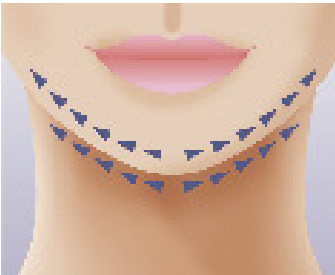Apply gel over the area to be worked.

Slide the handset using an upward circular motion from the centre outward, covering the entire area. Follow these simple instructions to help you get the best results.

Remember DO NOT apply SonixPlus on the neck and throat.

Max. Time	Energy	Type
3 min.	low-medium (1 to 7)	pulsed 

Chapter 4

For correct, safe use of SonixPlus

For correct and safe use of SonixPlus, you should keep the following application times in mind:

• **Per work area:** By "work area" we mean each part of the face where SonixPlus can be applied. The various work areas are:

- *Forehead:* **max. 3 minutes**
- *Eye contour:* **max. 3 minutes**
- *Cheeks and cheekbones:* **max. 3 minutes**
- *Lip contour:* **max. 3 minutes**
- ⊖ *Chin and jowls:* **max. 3 minutes**

• **Per daily session: 15 minutes maximum.**

By session we mean the application of SonixPlus on different work areas in the same day.

• **Between sessions: 48 hours minimum.**

By "between sessions" we mean the rest time between each session.

E.g.: In one week, a session schedule would be: Monday, Wednesday, Friday and Sunday.

• **Every 24 sessions (approx. 1 and a half months of treatment): rest 1 and a half months.**

After this rest period you can repeat the treatment.

VERY IMPORTANT

***** Never exceed the application times (a total of 15 minutes per session in one day).** Exceeding recommended applications can cause fatigue and skin irritation. On sensitive or allergic skins applications should be of short duration and at low energy settings.

***** And, above all, respect the rest period of 48 hours between sessions.** This is the time that your skin needs to rest.

Chapter 5

Frequently asked questions about SonixPlus

What is SonixPlus?

SonixPlus is a hi-tech device designed to stimulate skin tissue and improve the elasticity of the skin.

The technology used employs the emission of high frequency ultrasonic waves which help to improve the cell membrane metabolism, contributing to better micro-circulation. This technology increases the oxygenation of skin cells providing a completely safe and painless way to combat both the sagging skin and also remove wrinkles, dark circles, bags under the eyes, etc..

What are the advantages over other methods?

3MHz ultrasonic waves are safe, painless and non-invasive and they also aid sonophoresis, allowing the nutrients in skin care creams to penetrate deeper.

Is it painful?

It is completely painless and safe.

Is it normal to hear a slight buzzing in your ears while using SonixPlus?

Yes, some users may detect a slight buzzing sound during its application.

Can it be used on all parts of the body?

The 3Mhz SonixPlus technology has been specifically designed for use on the face. **Under no circumstances must SonixPlus be applied to the throat, neck, eyeballs, ears and ear canals.**

Is SonixPlus effective on "Crow's feet"?

Yes, SonixPlus is effective on "crow's feet", aiding the regeneration and tighten of wrinkle-causing cells.

Can SonixPlus be used to remove wrinkles?

Yes, SonixPlus is an effective way to remove wrinkles.

After how many sessions will I begin to see results?

You should notice results after the first five sessions, however results may vary depending on the skin of each individual.

How long should I wait between sessions?

It is very important to respect the 48 hour break between sessions, as this is the time that your skin needs to rest.

Do wrinkles reappear after using SonixPlus?

Although SonixPlus ensures effective and long-lasting results, wrinkles are part of the natural ageing process. To reduce and delay the onset of wrinkles the SonixPlus maintenance sessions are recommended.

Are there any side effects?

If you are healthy, SonixPlus has no reported side effects, nevertheless in some people it may produce a reddening of the skin that can last for a few minutes.

What is ultrasound gel for? When is it used?

The ultrasound gel is essential, ALWAYS use it before beginning a session. This will ensure that the ultrasonic waves penetrate correctly and prevent damage to the unit.

Can I use SonixPlus without applying gel?

Never, because the gel acts as a conductor for the ultrasonic waves on your face. If you do not use gel then the SonixPlus application loses its effectiveness and the machine can be severely damaged.

Can I use any type of gel?

We recommend using only TECNOVITA by BH gel, as the use of any other type may damage the machine.

Is it necessary to drink water?

No, it is not strictly necessary but it is advisable to hydrate yourself. The recommended quantity is at least 2 litres during treatment, as this helps to hydrate the skin.

When is it inadvisable to use SonixPlus?

- if you suffer from hyperpigmentation, broken capillaries
- if you have sunburn
- if you suffer from severe skin problems
- active skin infections
- facial paralysis
- If you are pregnant or think you might be
- If you are breast feeding.
- undergoing fertility treatments.
- couperosis
- temporal arteritis
- malignant tumours
- diabetes
- kidney problems
- physical or mental deficiencies
- metal implants or prostheses in the area of application
- scarring in the area of application
- physical malformations.
- epilepsy
- long-term treatment with medications: anti-inflammatories, anticoagulants, antibiotics or antihistamines.
- treatment for psychological or psychiatric disorders

Can I use SonixPlus if I am pregnant?

NEVER use SonixPlus if you are pregnant or suspect that you might be.

Can I use SonixPlus during lactation?

Do not use it if breastfeeding.

Can I use SonixPlus on the breasts?

Never apply it to the breasts.

Can I use SonixPlus if I have skin conditions?

No, It should not be used if you have any serious skin condition. Consult your doctor.

Can I use SonixPlus if I have ear implants?

No, it should not be used if you have hearing implants fitted.

Can I use SonixPlus if I have implants in my mouth, gold teeth, metal implants or similar?

Do not use it if you have these types of implant, as it may cause discomfort or tenderness in the area of application.

Can I use SonixPlus when I am bathing or showering?

Under no circumstances should the handset or the machine come into contact with water due to the risk of electrocution.

How do I clean SonixPlus and the handsets?

Never immerse SonixPlus or the handsets in water. Use a cloth slightly dampened with water.

Chapter 6

Important safety instructions

This product is supplied with a mains cable that has an earth wire and a plug with an earth pin. In the event that the appliance malfunctions then the earth connection will redirect the electrical current away from the machine and reduce the risk of an electrical shock.

DANGER! If the equipment is not earthed correctly then there is a risk of electrical shock.

Do not modify the electric plug supplied with the unit. If it does not fit the wall socket then get an electrician to install an electrical socket that does. This product must be used on a 220-240 volt circuit. Follow these recommendations:

1. Never use accessories that are not recommended by the manufacturer.
2. Do not use the machine if the mains cable becomes damaged or worn.
3. Keep the mains cable away from hot surfaces.
4. Do not use the machine outdoors.

Protection

This model is equipped with an electronic circuit breaker to protect both the electronic and electrical circuits.

This safety device (fuse) will be automatically triggered if the unit is subjected to abnormal force or heat.

Should this occur unplug the cable from the mains supply and change the fuse.

Warranty

We suggest that you read through this document carefully in order to fully understand the cover that the manufacturer offers as well as your responsibilities as owner.

GUARANTEES

GENERAL INFORMATION

EXERCYCLE, S.L. would like to thank you for purchasing this product and the trust you have placed in us. Our products have been subjected to rigorous quality checks. Nevertheless, should you experience any problem with this appliance during the guarantee period then please read through these warranty conditions and then get in touch with your local authorised dealer.

This guarantee certificate has been drafted in accordance with that stated in Directive 1999/44/EC of the European Parliament and of the Council 25th May, 1999, regarding on certain aspects of the sale of consumer goods and associated guarantees, incorporated into Spanish Law by virtue of Law 23/2003, 10th July; Guarantees on the sale of consumer goods.

To this effect, EXERCYCLE, S.L. VAT nº B-01029032, with business address c/ Zurrupitieta 22, Vitoria, Spain, guarantees this product under the following terms and conditions:

1. GUARANTEE PERIOD

The product covered by this guarantee certificate is protected for a period of TWO (2) YEARS as from the date of receiving the product, which will be taken as the date that appears on the sales receipt for the product or on the corresponding delivery docket should this be later.

It is essential that this be presented in the event of possible repairs. The policy regarding the warranty period may vary between different countries.

The stated warranty period is only applicable to Europe.

2. RIGHTS UNDER THE GUARANTEE

- Repair or replacement of the product which, in such instance, will be carried out within a reasonable period and without significant inconvenience to the consumer, taking into account of the nature of the goods and the purpose for which the consumer required them.
- Those cases in which the repair or replacement represents a disproportionate cost for the seller in comparison with other forms of remedy will be exempt assuming these do not prejudice the consumer.
- Repairs will be free for the consumer with regard to labour, materials and carriage.
- In the event that the machine is repaired or replaced with a new one, the guarantee period will be suspended temporarily until the new appliance has been delivered and will carry at least six months of guarantee, however under no circumstances will a new 2 year period be initiated.

The buyer may, as is apparent, require (i) a reduction in price or proceed to a (ii) recession of the contract when the repair or replacement is not carried out within a reasonable time or with significant inconvenience to the buyer. The reduction in price (i) will be proportional to the difference between the value that the goods would have had if they had been delivered in accordance with the contract and the actual value that the goods were worth at the time of the said delivery. Recession of the contract (ii) will not proceed when the non-conformity is of minor importance.

3. THE GUARANTEE DOES NOT COVER

- Defects attributable to use other than that established in the User Manual.
- Professional use or commercial leasing or hiring of the product.
- Poor installation, poor assembly or breakage of components due to improper installation, (the User Manual gives detailed instructions on how to install and use the machine).
- Improper use of the product (use without applying gel, exceeding operating times, exposure to high moisture levels, etc.) or problems arising from lack of maintenance (accumulation of dirt on any part of the product or its accessories, etc.). (The User Manual provides advice on how to install and use the product).
- Repairs carried out by persons not authorised by EXERCYCLE, S.L. or appliances that have been opened or tampered with by persons not authorised by the manufacturer.
- Deviations from the standard product, such as light vibrations, noises or breakage of components that do not affect the performance of the machine and which cannot be attributed to manufacturing faults.
- Scratches, scrapes, abrasions, knocks, falls, etc. caused by carelessness.
- Defects on the appliance caused by damage during transit.

4. CLAIMING AGAINST THE GUARANTEE


- Claims made for nonconformities will be limited to three years as from the date that the goods were delivered.

- The consumer must report the existence of the defect within a period of 2 months as from the date that the defect was discovered.
 - This guarantee does not affect the consumer's statutory rights under any applicable law, unless it is incompatible with the application of other actions derived from compensation sought for hidden faults in the contract of sale.
 - With regard to queries that may arise from the interpretation or application of this guarantee, these will be governed by the corresponding national law.
- Any discrepancy arising from the interpretation, validity or application of this contract will be subject to the courts and tribunals of Vitoria, Spain.

Certificates

By virtue of this document, Exercycle S.L. declares that this product meets the essential requirements and other relevant provisions of Directives 2004/108/EC and 2006/95/EC

Parts list

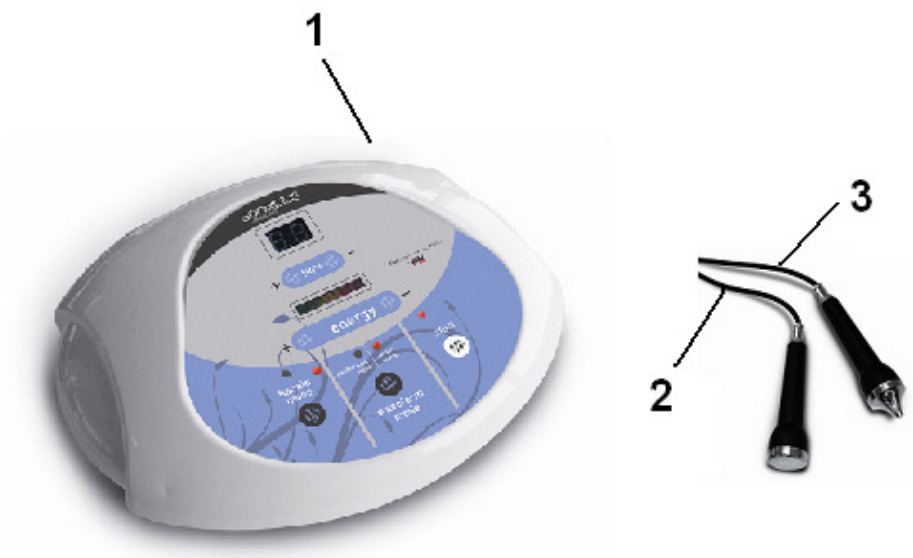- 1 - Main unit
- 2 - Flat handset
- 3 - Tapered handset

To order replacement parts, please indicate:

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BH RESERVES THE RIGHT TO MODIFY THE SPECIFICATIONS OF ITS PRODUCTS WITHOUT PRIOR NOTICE.

Do not hesitate to get touch with the Technical Assistance Service if you have any queries by phoning customer services (see last page in manual)

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